

## How can I get Roger?

**The answer is simple... you have to ask!**

Reasonable educational accommodations are granted upon student request. Your rights are protected by Section 504 of the Rehabilitation Act of 1973 and the Americans with Disabilities Act of 1990.

Inquire with your university's disability coordinator or accessibility center to learn how you can acquire products like Phonak Roger.

Please visit <http://us.morethanahearingaid.com/college/> for more information.

Audiologist: \_\_\_\_\_

Contact # / Email: \_\_\_\_\_

Current Roger System

Mic(s): \_\_\_\_\_

Receiver: \_\_\_\_\_

Recommended Roger System

Mic(s): \_\_\_\_\_

Receiver: \_\_\_\_\_

## The Roger Difference

**What is Roger?**

Roger is a digital wireless standard that uses 2.4 GHz tech to boost speech understanding in all noise levels and over distance. It wirelessly transmits the speaker's voice directly to the listener, offering speech-in-noise improvements of up to 35% over previous Phonak wireless technology and 54% over other FM systems.

**Features that set Roger above any other system on the market:**

1. Ease of use – no more channel planning, this is done automatically – if interference is detected in the environment the Roger system will select a new channel automatically.
2. Eavesdrop security – encrypted signal.
3. MultiTalker Network – ability to pair together multiple mics to cover a larger room or multiple speakers!
4. Continuously monitors background noise and speech in environment – noise suppression and gain settings adjust automatically to provide the best possible speech understanding.
5. Bluetooth capable – connect to computers, cell phones etc.
6. Fast transmission time – minimal time delay, no echo.
7. Compatible with virtually every hearing aid/cochlear implant; can also be utilized with headphones or Roger Focus.
8. Full audio frequency bandwidth to ensure speech understanding across the speech spectrum.
9. MultiBeam Technology - by utilizing multiple microphones in six directions, speech from a 360° radius is calculated and compared. The direction with the best clarity is automatically selected.

## Life is on

At Phonak, we believe that hearing well is essential to living life to the fullest. For more than 70 years, we have remained true to our mission by developing pioneering hearing solutions that change people's lives to thrive socially and emotionally. Life is on.

<http://us.morethanahearingaid.com/college/>



## Assistive Listening Devices for college

Focus on class rather than on hearing

**sonova**  
HEAR THE WORLD

028-1919-03/2018-07 /J5/Printed in US © 2018 Sonova AG All rights reserved

A Sonova brand

**PHONAK**  
life is on

# When you can focus on class rather than on hearing, life is on

Even the most powerful hearing aids have limitations. The Roger portfolio consists of wireless microphones and receivers. They pick up the speaker's voice and transmit it to the receivers connected to your hearing aids, so you can hear speech clearly over distance and in noise.

## Roger microphones



**Roger Table Mic II**  
Designed for various discussions; several microphones can be connected to cover large groups.



**Roger Select™**  
A versatile microphone which automatically selects the person who is talking and with Bluetooth® for phone calls.



**Roger Pen™**  
A handy microphone with portable design and Bluetooth for phone calls.

## Roger receivers



**Roger design-integrated**  
Click onto Phonak hearing aids and selected implants.



**Roger X**  
Universal miniature receiver for behind-the-ear aids, implants and streamers.



**Roger MyLink**  
Neck loop receiver for hearing aids or implants with a telecoil.

Situations with background noise or where the sound source is at a distance are ones where Roger truly excels:

### Small discussions




Small discussions can be overwhelming for students wearing hearing aids. Place one microphone on the table and hear your classmates clearly.

### Large discussions



In a large group discussion, you might need more than one table microphone. You can use as many Roger Table Mic II as necessary.


### Lecture with a professor

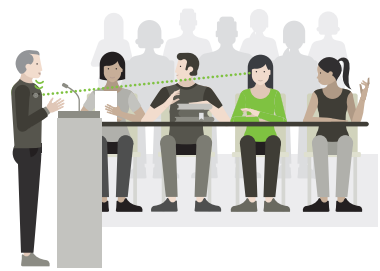

Professors often stand at a distance from the audience. When the presenter wears Roger Select, it works smoothly together with the table microphone.


Additional situations where Roger can help overcome communication challenges:

### One-on-one conversation in noise

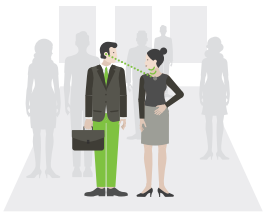


When surrounded by distracting noise hang a Roger microphone around your classmate's neck or clip it on to the collar.

### Multimedia

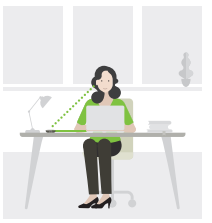

Video conferences and webinars are frequently used in distance learning. Simply connect Roger Table Mic II, using the audio cable.

### Phone calls

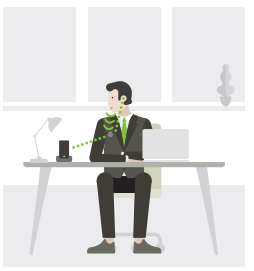

By connecting Roger Pen or Roger Select to your phone via Bluetooth you can hear phone calls directly through your hearing aids or implants.

## Bridging the understanding gap

### What is Roger?

Roger is a digital wireless standard specifically designed to help boost your hearing performance<sup>1</sup>.

### Communicate

The modern classroom can be a challenging listening environment. It is also a place where successful communication is key in order to succeed.

### Participate

When affected by hearing loss, you become a master at finding ways around difficult hearing situations. This can be time consuming, emotionally demanding and often limits your ability to fully participate and engage with your classmates.

### Contribute

Roger wireless microphones can be used together with modern hearing aids, a solution that boosts your performance and makes it easier for you to interact in the classroom.

<sup>1</sup>Mülder, H. (2013). Roger and hearing instruments. Excellent speech understanding in high noise levels. Phonak Field Study News, retrieved from [www.phonakpro.com/evidence](http://www.phonakpro.com/evidence), accessed November 14th, 2017.